


ITALIANO (Lingua originale)	AUTOCLAVE DI STERILIZZAZIONE A VAPORE D'ACQUA
ENGLISH (translation from the original in Italian)	WATER STEAM STERILIZATION AUTOCLAVE
FRANÇAIS (traduction de l'original en italien)	AUTOCLAVE DE STÉRILISATION À LA VAPEUR D'EAU
DEUTSCH (Übersetzung aus dem Original in Italienisch)	WASSERDAMPF-STERILISATOR
ESPAÑOL (traducción del original en italiano)	AUTOCLAVE DE ESTERILIZACIÓN POR VAPOR DE AGUA
POLSKI (przekład z oryginału w języku włoskim)	AUTOKLAW DO STERYLIZACJI PAROWEJ
TÜRKÇE (İtalyan orijinalinden çeviri)	SU BUHARLI STERİLİZASYON OTOKLAVI
NORSK (oversettelse fra originalen på italiensk)	AUTOKLAV FOR STERILISERING MED VANNDAMP
ČESKY(překlad z originálu v italštině)	AUTOKLÁV PRO STERILIZACI VODNÍ PAROU
NEDERLANDS (vertaling uit het Italiaans)	AUTOCLAAF VOOR STERILISATIE MET WATERDAMP
ELLHNIKA (αση από το πρωτότυπο στα ιταλικά)	ΑΥΤΟΚΑΥΣΤΟΣ ΚΛΙΒΑΝΟΣ ΑΠΟΣΤΕΙΡΩΣΗΣ ΜΕ ΥΔΡΑΤΜΟΥΣ
РУССКИЙ (Перевод с оригинала в итальянском)	СТЕРИЛИЗУЮЩИЙ АВТОКЛАВ НА ВОДЯНОМ ПАРУ
PORTUGUÊS (tradução do original em italiano)	AUTOCLAVE DE ESTERILIZAÇÃO A VAPOR D'ÁGUA
SVENSKA (översättning från originalspråket italienska)	AUTOKLAV FÖR STERILISERING MED VATTENÅNGA
SUOMI (käännös alkuperäisestä italialaisesta teksistä)	VESIHÖYRYLLÄ TOIMIVA STERILOINTIAUTOKLAAVI
اللغة العَرَبِيَّة (ترجمة من اللغة الإيطالية)	*جهاز أوتوكلاف للتعقيم ببخار الماء*
MAGYAR (az eredeti fordítása)	VÍZGŐZÖS STERILIZÁLÓ AUTOKLÁV
БЪЛГАРСКИ (превод от италиански)	АВТОКЛАВ ЗА СТЕРИЛИЗАЦИЯ С ВОДНА ПАРА
DANSK (Oversat fra italiensk)	AUTOKLAVE TIL VANDDAMPSTERILISERING
LIETUVIŲ K. (versta iš italų k.)	STERILIZAVIMO VANDENS GARAIS AUTOKLAVAS
ROMÂNĂ (traducere din limba italiană)	AUTOCLAVĂ DE STERILIZARE CU VAPORI DE APĂ
SLOVENSKY (preklad talianského originálu)	AUTOKLÁV SO STERILIZÁCIOU VODNOU PAROU
LATVISKI (tulkojums no angļu valodas)	ŪDENS TVAIKA STERILIZĀCIJAS AUTOKLĀVS
意大利语（原始语言）	**水蒸汽**杀菌机
УКРАЇНСЬКА (переклад з англійської)	ВОДНО-ПАРОВИЙ СТЕРИЛІЗАЦІЙНИЙ АВТОКЛАВ
EESTI KEEL (tõlgitud inglise keelest, originaalkeel on itaalia keel)	VEEAURU STERILISEERIMISAUTOKLAAV
SRPSKI (prevod iz originala na talijanski)	AUTOKLAV ZA STERILIZACIJU NA VODENU PARU

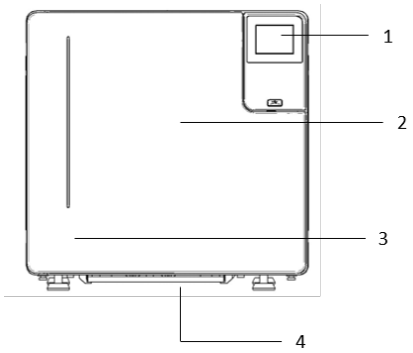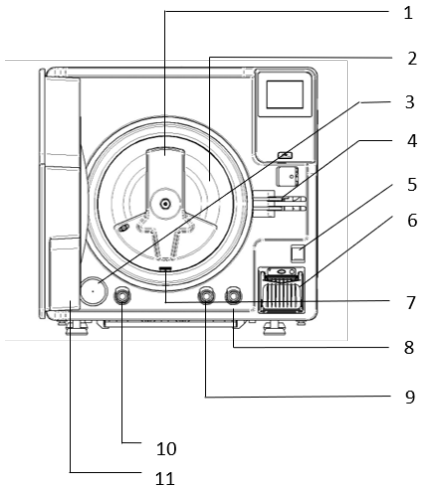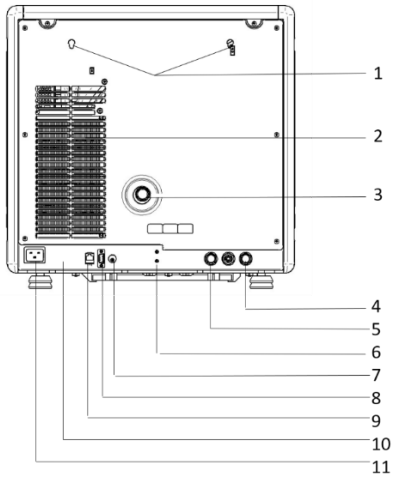NOTA

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NOTE

COMPLETE USER MANUAL AVAILABLE ON SUPPLIED USB DEVICE.

Ref.	Description
1.	TOUCH control panel
2.	Door
3.	Model Label
4.	Dust filter

Ref.	Description
1.	Steam diffuser
2.	Sterilization chamber
3.	Bacteriological filter
4.	Door locking system
5.	Power switch
6.	Printer (optional)
7.	Chamber drain filter
8.	Distilled water top-up quick coupling
9.	Distilled water drainage quick coupling
10.	Waste water drainage quick coupling
11.	Door

Ref.	Description
1.	Fastening slots for rear spacers
2.	Heat exchanger
3.	Safety valve
4.	Connection for direct water drainage
5.	Connection for automatic distilled water filling (only for PURE 100 / 500, EV AUX kit and external pump kit)
6.	Network fuse fastening screws
7.	Automatic filling electrical connection (only for PURE 100 / 500, EV AUX kit and external pump kit)
8.	Serial cable connection
9.	Ethernet cable connection (max. length 29 m)
10.	Data plate / SERIAL NUMBER LABEL
11.	Power cable connection

CYCLES	TYPE OF LOAD	MAX KG				
		S17	S22	17	22	28
B 134° C Universal B 121° C Universal B 134° C Prion (*)	Unwrapped porous materials	-	-	1,00	1,20	1,50
	Porous materials in single package	-	-	0,75	1,00	1,25
	Porous materials in double package	-	-	0,60	0,75	0,90
	Solid and hollow materials in single package	-	-	3,00	4,00	5,00
	Solid materials and cables not packed	-	-	6,00	7,50	9,00
	Solid and hollow instruments in double package	-	-	1,50	2,00	2,50
S 134° C Veloce	Unwrapped solid and hollow instruments	-	-	2,00	2,00	2,00
S 134° C Solid	Unwrapped solid and hollow instruments “B”	6,00	7,50	6,00	7,50	9,00
	Solid and hollow instruments “B” in single package	3,00	4,00	3,00	4,00	5,00
S 121° C Solid S 134 ° C Prion (*)	Solid and hollow instruments “B” in single package	3,00	4,00	-	-	-
N 134° C Solid N 121° C Solid	Unwrapped solid and hollow instruments “B”	6,00	7,50	-	-	-

(*) NOT AVAILABLE FOR US MARKET

NOTE

Definition of hollow loads in accordance with standard EN13060:2023

In this manual the term "hollow loads" refers both to the elements defined as "narrow lumen" (paragraph 3.18 EN 13060:2023) and to the elements defined as "simple hollow" (paragraph 3.30 EN 13060:2023).

The term "hollow load B" refers ONLY to the elements defined as "simple hollow" (paragraph 3.30 EN 13060:2023).

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